


MICROSURGERY KNIVES

MIKROCHIRURGISCHES EINMALMESSER
BISTOURI POUR MICROCHIRURGIE JETABLE
BISTURI MONOUSO PER MICROCHIRURGIA
CUCHILLETE DE MICROCIURUGIA DESECHABLE
BISTURI DE MICROCIURURGIA DESCARTÁVEL
TEK KULLANIMLIK MİKROCERRAHİ BİÇAK



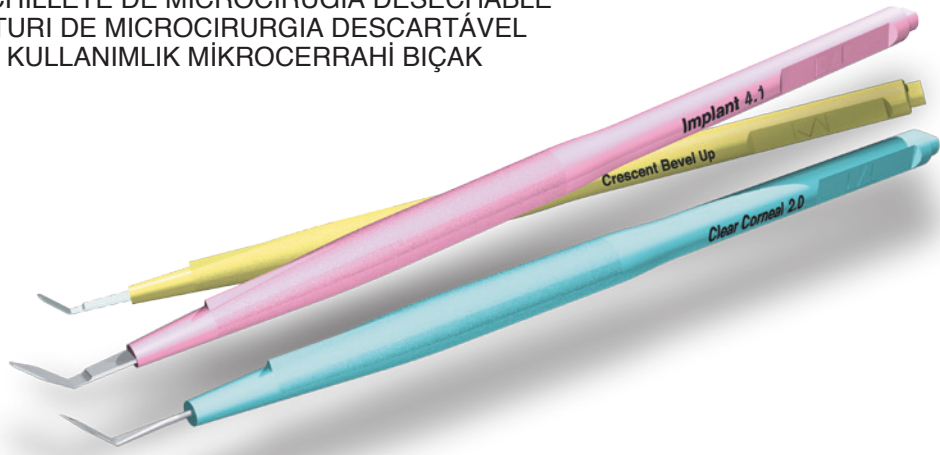


STAB KNIVES -----	3	page
CRESCENT KNIVES -----	4	
BLEB KNIVES -----	4	
SLIT KNIVES -----	5	
CLEAR CORNEAL KNIVES -----	6	
TRAPEZOID KNIVES -----	7	
IMPLANT KNIVES -----	7	
MVR KNIVES -----	8	
PARACENTESIS KNIVES -----	9	
SIDEPORT KNIVES -----	10	

Regular Type



Sterile and ready to use
Each knife in a protection tray
5 pieces per dispenser box
GS1-128 bar-code

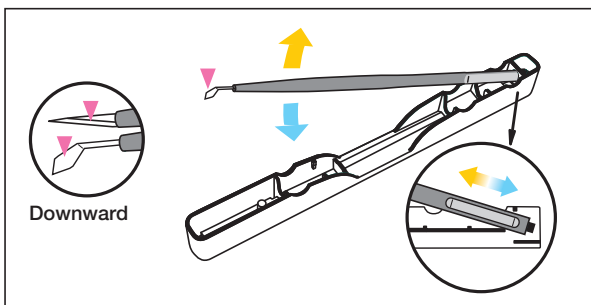


Safety Type



Sterile and ready to use
Blade covered with protective sheath
5 pieces per dispenser box
GS1-128 bar-code


Blade protection allows safety handling and avoids accidental injuries
Precisely engineered design for easy twist and glide action with both hands

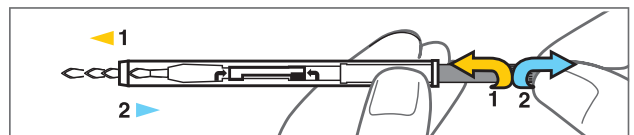

STAB KNIVES

Cataract Surgery / Scleral Tunnel Incision - Pregroove, Side port
/ Clear Corneal Incision - Side port


Precisely ground in 3 bevels for ultra sharp incision and piercing.

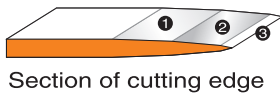



Section of cutting edge

Distinction of Kai's Stab Knife

Stab Knife blade is made of razor blade with 3 bevel edge-setting.
The razor edge-setting method ensures the consistent sharpness and precision of Stab Knife.

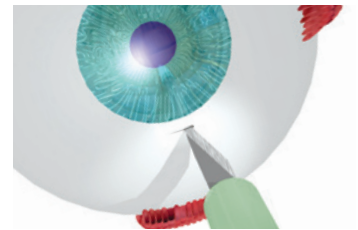

Regular Type

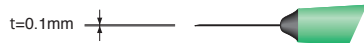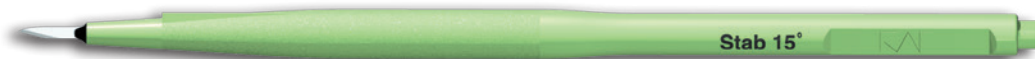

A-0F angle 0°		A-15F w=1.7mm angle 15°		A-225F w=1.7mm angle 22.5°	
A-30F w=1.7mm angle 30°		A-45F w=1.7mm angle 45°			

Safety Type




SK-A15 w=1.7mm angle 15°		SK-A225 w=1.7mm angle 22.5°		SK-A30 w=1.7mm angle 30°	
---------------------------------------	--	--	--	---------------------------------------	--

CRESCENT KNIVES

Cataract Surgery / Scleral Tunnel Incision - Scleral Tunnel

Either straight or tapered sides are available to create the optimal tunnel or flap.

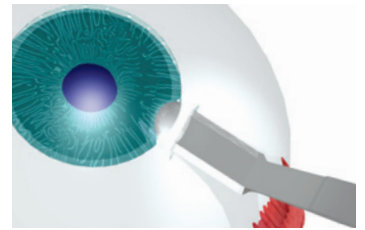

Regular Type



Anti-glare type Bevel Up

CPU-16AGF w=1.6mm angle 45°	CPU-22AGF w=2.2mm angle 45°	CPU-26AGF w=2.6mm angle 45°
CWU-26AGF w=2.6mm angle 60°		

Safety Type



Anti-glare type Bevel Up

SK-CPU16 w=1.6mm angle 45°	SK-CPU22 w=2.2mm angle 45°	SK-CPU26 w=2.6mm angle 45°
---	---	---

BLEB KNIVES

Glaucoma / Trabeculectomy / Needling

Regular Type


Anti-glare type

BKS-10AGF w=1.0mm Straight	BKB-10AGF w=1.0mm Bent Bevel Up angle 45°	BKB2-10AGF w=1.0mm Bent Bevel Up angle 45°
---	--	---

SLIT KNIVES

Cataract Surgery / Scleral Tunnel Incision - Anterior Chamber Entry

Many sizes are available to ensure the exact incision for your choice of phaco tip and operating access. Precision created by computer driven grinding technology.

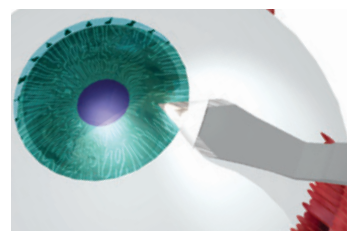

Regular Type



Anti-glare type

Bevel Up


SLU-13AGF w=1.3mm	SLU-14AGF w=1.4mm	SLU-15AGF w=1.5mm
SLU-16AGF w=1.6mm	SLU-17AGF w=1.7mm	SLU-18AGF w=1.8mm
SLU-19AGF w=1.9mm	SLU-20AGF w=2.0mm	SLU-22AGF w=2.2mm
SLU-23AGF w=2.3mm	SLU-24AGF w=2.4mm	SLU-25AGF w=2.5mm
SLU-265AGF w=2.65mm	SLU-275AGF w=2.75mm	SLU-28AGF w=2.8mm
SLU-285AGF w=2.85mm	SLU-30AGF w=3.0mm	SLU-32AGF w=3.2mm

Safety Type



Anti-glare type

Bevel Up

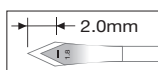


A knife with "M" at the end of Ref number has a line and maximal width mark at 2.0 mm from the blade tip.

SK-SLU18M w=1.8mm	SK-SLU20M w=2.0mm	SK-SLU22 w=2.2mm
SK-SLU23 w=2.3mm	SK-SLU24 w=2.4mm	SK-SLU28 w=2.8mm
SK-SLU30 w=3.0mm	SK-SLU32 w=3.2mm	

CLEAR CORNEAL KNIVES

Cataract Surgery / Clear Corneal Incision - Whole Incision

The latest technology, chemical etching process creates this ultra sharp blade, producing a double bevel perfectly suited to clear corneal incisions. The new geometry allows high maneuverability for crafting the phaco self-sealing corneal tunnel.

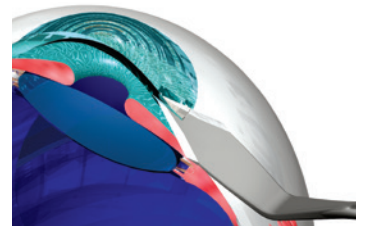

Regular Type

Anti-glare type

Double Bevel

CCR-18AGF w=1.8mm 	CCR-20AGF w=2.0mm 	CCR-22AGF w=2.2mm 
CCR-23AGF w=2.3mm 	CCR-24AGF w=2.4mm 	CCR-25AGF w=2.5mm 
CCR-265AGF w=2.65mm 	CCR-275AGF w=2.75mm 	CCR-28AGF w=2.8mm 
CCR-30AGF w=3.0mm 	CCR-32AGF w=3.2mm 	

Safety Type

A knife with "M" at the end of Ref number has a line and maximal width mark at 2.0 mm from the blade tip.

Anti-glare type

Double Bevel

SK-CCR18M

SK-CCR20M

SK-CCR22

SK-CCR23

SK-CCR24

SK-CCR28


SK-CCR30

SK-CCR32


TRAPEZOID KNIVES

Cataract Surgery / Clear Corneal Incision - Whole Incision

The Trapezoidal shape creates a tunnel with graduated sides, non parallel, to provide more maneuverability for the phaco tip.

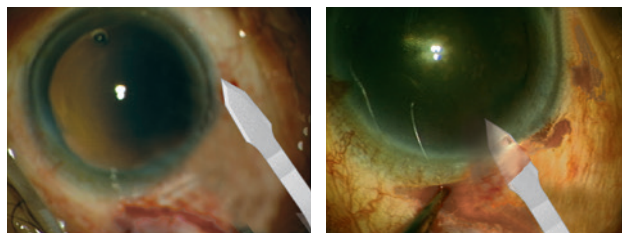


Regular Type

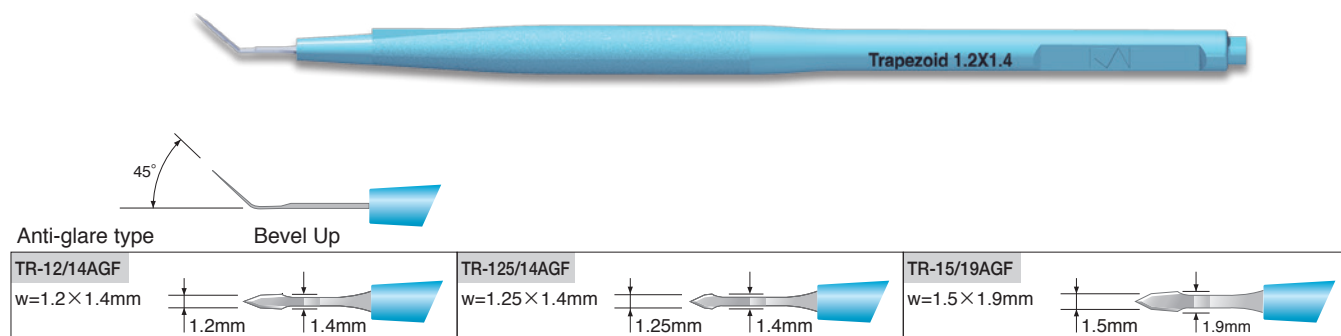


IMPLANT KNIVES

Cataract Surgery / Scleral Tunnel Incision - Incision Widening for IOL

Manufactured to extend the phaco tip incision to create the exact size required to implant your intra ocular lens of choice.

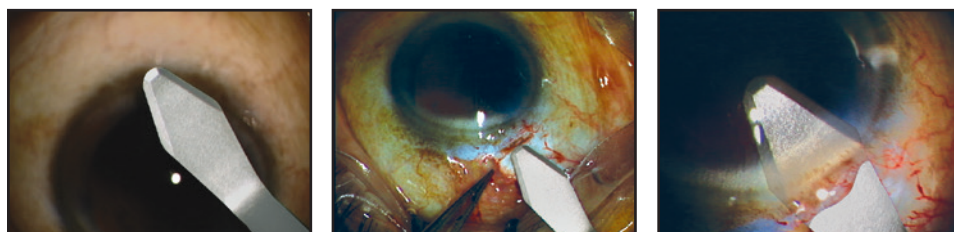


Regular Type

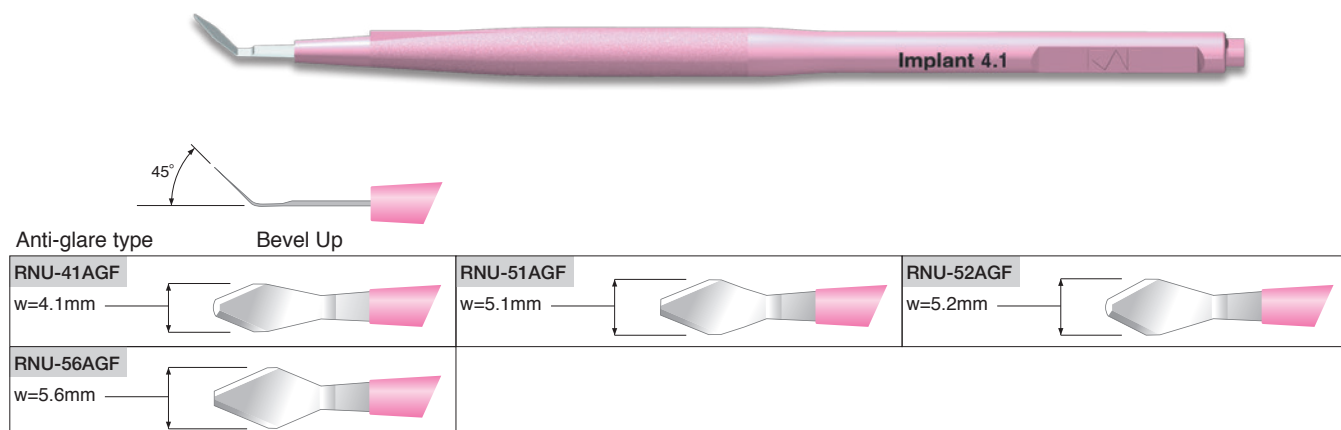

MVR KNIVES

Cataract Surgery / Scleral Tunnel Incision - Side port
/ Clear Corneal Incision - Side port
/ Retina Surgery - Scleral Incision

Specifically designed for very precise vitrectomy surgery.
Available in both straight and angled styles.

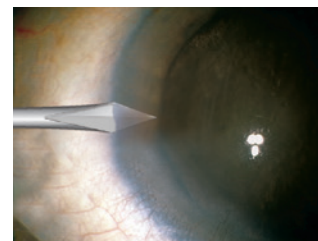

Regular Type



VK-25GF $\Phi=0.5\text{mm}$ $\Phi=0.5\text{mm}$ Straight	VK-23GF $\Phi=0.6\text{mm}$ $\Phi=0.6\text{mm}$ Straight	VK-20GF $\Phi=0.9\text{mm}$ $\Phi=0.9\text{mm}$ Straight
VK-19GF $\Phi=1.0\text{mm}$ $\Phi=1.0\text{mm}$ Straight		
VKB-20GF $\Phi=0.9\text{mm}$ $\Phi=0.9\text{mm}$ Bent		

Safety Type




SK-VKB20G $\Phi=0.9\text{mm}$ $\Phi=0.9\text{mm}$ Bent
--

PARACENTESIS KNIVES

Cataract Surgery / Scleral Tunnel Incision - Side port
/ Clear Corneal Incision - Side port

Produced by grinding technology for highly accurate second instrument access. Post operatively, completely self sealing.

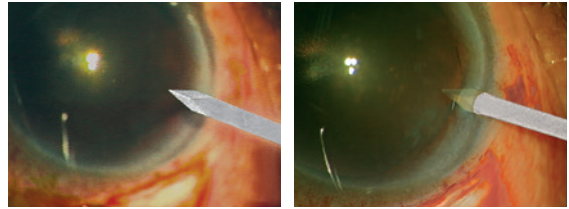

Regular Type



Anti-glare type

PCT-10F

w=1.0mm



Straight



PCB-10AGF

w=1.0mm



angle 45°



SIDEPORT KNIVES




Cataract Surgery / Scleral Tunnel Incision - Side port
/ Clear Corneal Incision - Side port

Regular Type

SIDEPORT KNIFE TRAPEZOID



Suitable for Eyes with Shallow Anterior Chambers

Compared to the 15 degree or 22.5 degree stab knives, Sideport Knife Trapezoid can create trapezoidal incision with short incision distance. It minimizes the risk of damaging nucleus and especially effective to the patients with shallow anterior chambers. Also, the incision size can be adjusted by the depth of penetration, and it is easy to enlarge the incision with the sufficient length of the blades on both sides.

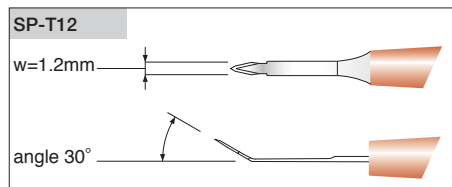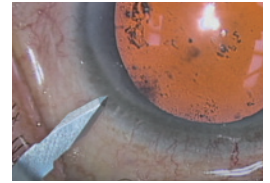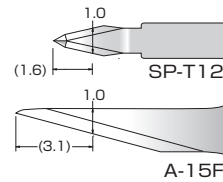

Improved Maneuverability


It is angled at 30 degree, and the head is compact. These features make penetration into the eyes at the ideal angle possible even for the patients with deep set eyes or narrow palpebral fissures.

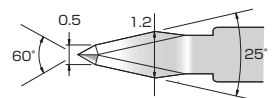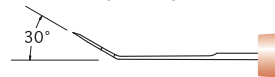


Regular Type

SIDEPORT KNIFE STRAIGHT

A 1 mm-width sideport can be stably prepared. The 15°-angle-blade allows easy expansion of incision.

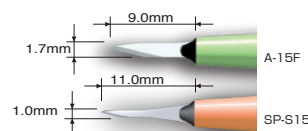




It can be used in the same way as a stab knife.

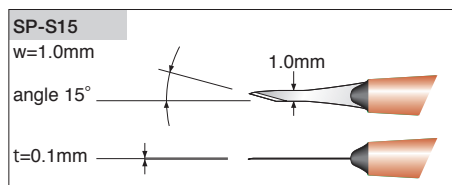

Regular Type

SIDEPORT KNIFE UTILITY

Able to make a stable sideport.
Can easily create a port when using a retractor.
Able to incise lamellar sclera with a single sideport knife utility.

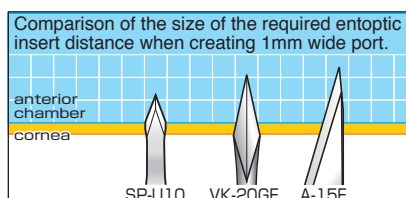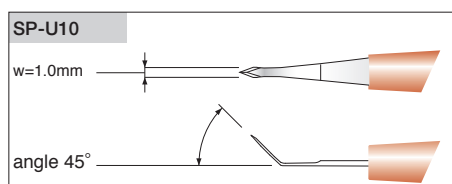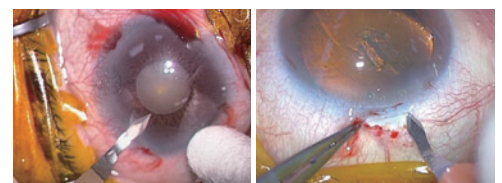

MICROSURGERY KNIVES

MIKROCHIRURGISCHES EINMALMESSER
BISTOURI POUR MICROCHIRURGIE JETABLE
BISTURI MONOUSO PER MICROCHIRURGIA
CUCHILLETE DE MICROCHIRURGIA DESECHABLE
BISTURI DE MICROCHIRURGIA DESCARTÁVEL
TEK KULLANIMLIK MIKROCERRAHI BİÇAK



kai industries co., ltd.

1110, Oyana, Seki City, Gifu Pref., 501-3992, Japan
Phone: +81 (0) 575 28 6600
Facsimile: +81 (0) 575 28 6611
E-mail: info-ind@kai-group.com
www.kai-ind.co.jp/en

CE
0197

Product specifications may change without notice.



kai Europe GmbH

Kottendorfer Straße 5, 42697 Solingen, Germany
Phone: +49 (0) 212 23 23 80
Facsimile: +49 (0) 212 23 23 899
www.kai-europe.com

KIM-021 0608/1407 K 5K9