


677PMY

Single piece, sterile, foldable acrylic aspheric hydrophilic posterior chamber intraocular lens (IOL) with UV blocker and blue light filter. The optic is trifocal. The IOL is preloaded, which is part of an injection system, delivered in a special lens holder, that must be inserted into a dedicated injector before implantation.

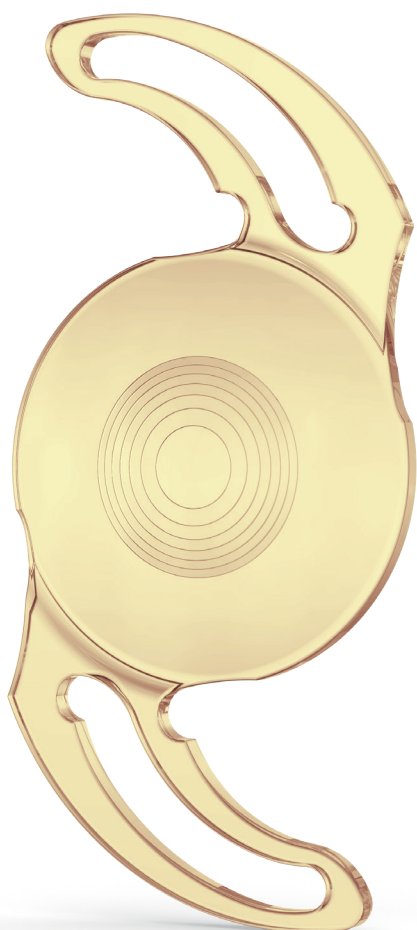

YELLOW IOL

POWER RANGE

Standard Powers

+8.0 D → +30.0 D (0.5 D steps)*

Addition

+3.5 D near
+1.75 D intermediate

MATERIAL

Copolymer of hydrophobic and hydrophilic monomers with 25% water content, UV blocker and blue light filter

Refractive Index

1.46

Abbe Number

58

DESIGN

Optic Diameter
Overall Diameter

6 mm

13 mm

Optic Design

Aspheric - Aberration Neutral
Biconvex (+8.0 D → +35.0 D)

Nominal A-constant

118.9

Diffractive Zone

3 mm diameter diffractive array on the
anterior surface utilizing EPS technology

Haptic Design

Posterior vaulting fenestrated C-loops
with 0° angulation

PCO Prevention

360° sharp edge

Shelf Life

3 years after sterilization

INJECTION

Estimated Incision Size

2.2 - 2.4 mm

Compatible Injector

MEDJET PIL-MA SINGLE-USE INJECTOR

CONSTANTS

for optical and immersion ultrasound measurements

SRK/T A-Constant

A = 118.828**
(optimized)

Haigis

a0 = 0.19**
a1 = 0.192**
a2 = 0.173**

Hoffer Q

pACD = 5.431**

Holladay I

SF = 1.682**

Barrett Universal II

LF = 1.79

Holladay II

ACD = 5.45
SF = 1.69


*Extreme powers available: +30.5 D → +35.0 D (0.5 D steps)

**Optimized constants based on more than 400 patients

677PMY

PRELOADED


Trifocal
LIBERTY[®]
MedJet PIL-MA

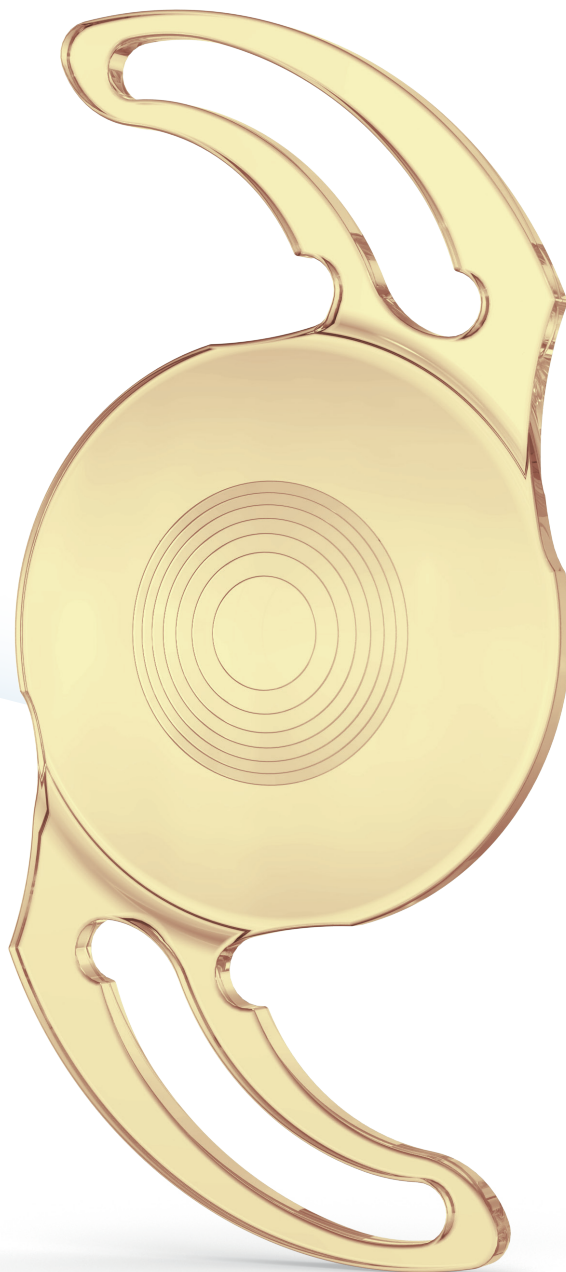

TRIFOCAL



PRELOADED

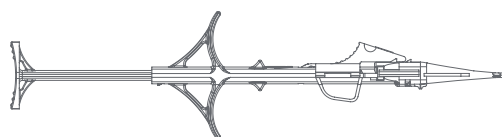

MEDICONTUR

Medicontur

Medical Engineering Ltd
export@medicontur.com
www.medicontur.com

Head Office

Herceghalmi út 1.
2072 Zsámbék
Hungary

International Office

Chemin du Champ-des-Filles 36
1228 Plan-les-Ouates
Switzerland / Geneva